


Measuring & Sizing Guide



1 Band Size

Measure the circumference just under the breast for band size

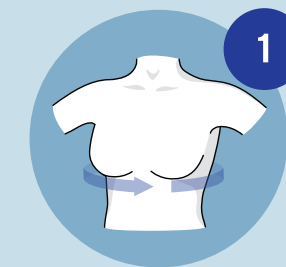


2 Cup Size

Measure the circumference of the fullest part of the breast for cup size

Sizing Chart		XS	S	M	L	XL
① Band Size	CM	68 - 77 cm	78 - 87 cm	88 - 97 cm	98 - 107 cm	108 - 117 cm
	IN	26.8 - 30.3 in	30.7 - 34.3 in	34.6 - 38.2 in	38.6 - 42.1 in	42.5 - 46 in
② Cup Size	AA/A	80 - 89 cm 31.5 - 35 in	90 - 99 cm 35.4 - 39 in	100 - 109 cm 39.4 - 43 in	112 - 119 cm 44.1 - 46.9 in	120 - 129 cm 47.2 - 50.8 in
	PIP CODE	4465316	4465340	4465332	4465324	4465308
	B/C	84 - 93 cm 33.1 - 36.6 in	94 - 103 cm 37 - 40.6 in	104 - 113 cm 41 - 44.5 in	116 - 123 cm 31.5 - 45.7 in	124 - 133 cm 48.8 - 52.4 in
	PIP CODE	4465225	4465233	4465241	4465209	4465217
	D/E	88 - 97 cm 34.6 - 38.2 in	98 - 107 cm 38.6 - 42.1 in	108 - 117 cm 42.5 - 46.1 in	120 - 127 cm 47.2 - 50 in	128 - 137 cm 50.4 - 54 in
	PIP CODE	4465274	4465290	4465282	4465258	4465266

*Rounding may apply to inch measurements

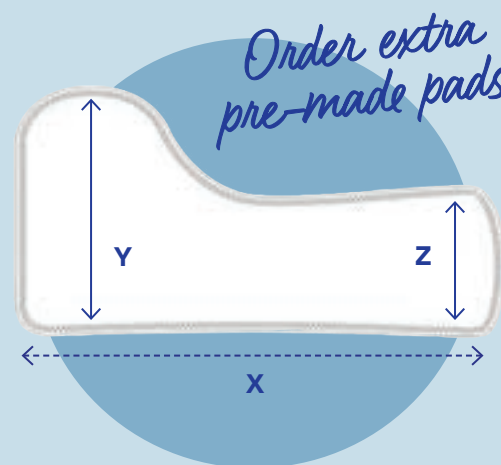


Item Number	Size	X	Y	Z
T57135120199970	XS	260 mm	155 mm	100 mm
T57135120299970	S	285 mm	160 mm	100 mm
T57135120399970	M	310 mm	165 mm	105 mm
T57135120499970	L	335 mm	170 mm	110 mm
T57135120599970	XL	360 mm	180 mm	120 mm

The pre-made Mobiderm Breast Pads accommodate all cup sizes and come in Small Block only.

Available from thuasne-store.co.uk

Or you can ...

Order additional **Mobiderm® Mobilising Pads** available in Big and Small Block (sold separately or available on prescription) to cut and use for different shapes or sizes

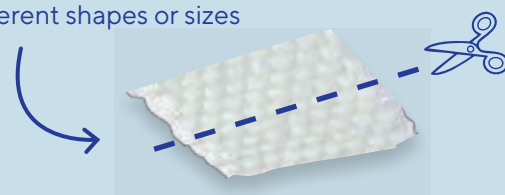

Pad Size	Block Size	PIP CODE
25cm x 1m	Small Block 5mm x 5mm	4197794
20cm x 1m	Big Block 15mm x 15mm	4197786
25cm x 25cm	Small Block 5mm x 5mm	4197802



Care Instructions

MOBIDERM® INTIMATE BRA AND PADS:

- Gentle machine wash on cold setting.
- Use mild detergents. Avoid use of bleach, fabric softener or other harsh chemicals.
- Place the garment in a wash bag with similar colors.
- Hang or lay flat to dry the bra, and keep it away from direct sunlight.
- Ensure all clasps and fasteners are secure prior to washing.
- Do not iron.



Materials

PRODUCT CONTAINS:

- 64% Nylon
- 24% Elastane
- 7% Polyester
- 5% Polyurethane

Soft & Breathable



Indications

- Acute or chronic oedema and/or lymphoedema of the breast and/or breast wall
- Postoperative oedema or oedema after radiotherapy (only on intact skin)
- Chronic infiltration after reconstruction, remaining chronic breast wall oedema



Black version
COMING
SOON
in 2027



Scan the QR code
to learn more about the
Mobiderm® Intimate Bra

Thuasne UK Ltd
Unit 4 Orchard Business Centre
North Farm Road
Tunbridge Wells, TN2 3XF

Phone 01892 481620
Email orders@thuasne.co.uk
thuasne.co.uk
thuasne-store.co.uk



THUASNE®

MOBIDERM® Intimate Bra

For Breast and Chest
Wall Oedema

**NOW AVAILABLE
ON PRESCRIPTION**



MOBIDERM® Intimate Bra

The Mobiderm® Intimate Bra provides optimum care for breast related oedema or lymphoedema by providing coverage on the breast, thoracic wall and back.


The Mobiderm® Intimate Bra provides strong compression and integrates removable and discreet foam compression pads to create high and low pressure differentials on the affected regions, in result, softening tissue and improving lymph drainage. Moreover, this controlled compression solution helps minimise scar tissue and the development of hematomas. The Mobiderm Intimate Bra is the bra you'll want to wear all day!

-  **Manufacturing Type:** Off-The-Shelf
-  **Time of Day:** Day and Night
-  **Compression Level:** 5-10mmHg
-  **Color:** Nude



The MOBIDERM® Technology


Mobiderm® foam blocks provide a comfortable, massaging effect to support clinical efficacy and compliance. This technology is made up of big or small foam blocks encased between two non-woven bandages. Due to the difference in pressure created between the support area of the blocks and their surrounding area, it promotes the flow of lymphatic fluid and optimises the drainage efficacy. ⁽²⁾

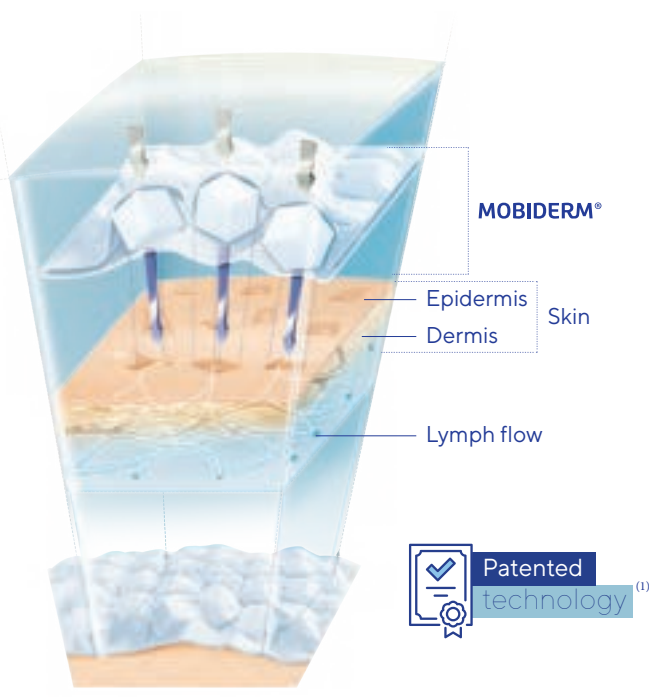

What to Love About the Mobiderm® Intimate Bra?

Proven Efficacy of Mobiderm® Technology

- ✓ Reinforce the effects of the daytime treatment ⁽²⁾⁽³⁾
- ✓ Designed to maintain the volume reduction of the oedema ⁽³⁾ and lymphoedema
- ✓ Promotes the flow of lymphatic fluid and optimises the drainage efficacy ⁽²⁾
- ✓ Mobilises oedematous or subcutaneous tissue (e.g. fibrosis)

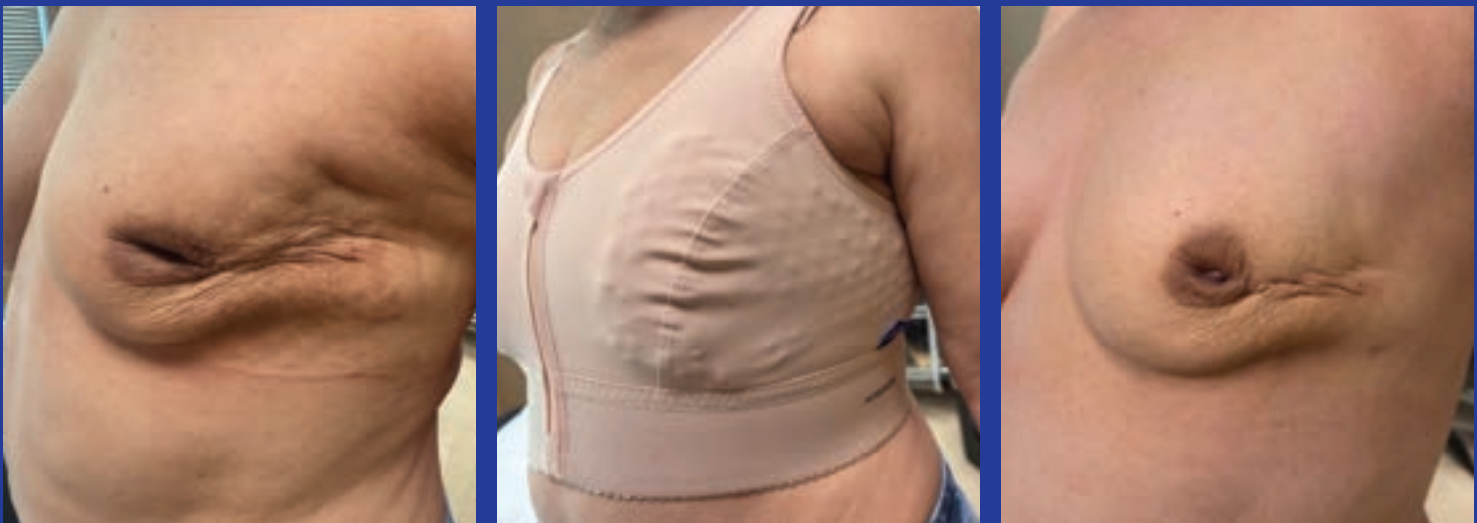

THERAPIST TESTIMONIAL | Linda Roherty, PT, CLT (USA)
"Oh my gosh! My patient is so very happy! Her breast has improved significantly, and she is no longer experiencing breast pain or nipple burning. Thank you!"

Mobiderm Intimate Bra



MOBIDERM BREAST PADS

Included in the box!



- ✓ 2 Small Block Breast Pads included
- ✓ Fits into the integrated pockets
- ✓ Provides discreet low profile therapy

Spare pads also available from thuasne-store.co.uk